ALPINGLOW SUMMIT TRAVELER

User Manual

Enjoy the glowing wellness, anywhere you go!

Introduction

- 1.SUMMIT TRAVELER is a portable and rechargeable led therapy light, with built-in battery, you can enjoy your light therapy anywhere and anytime. Slip the this SUMMIT TRAVELER into your backpack, bring it to the office or gym. Take it on the car, or on the plane.
- 2.SUMMIT TRAVELER adopt 15pcs 5w double chip LEDs and 30 degree lens. it is biggest high-power led chip, greatly improved the irradiance and light efficiency.
- 3.SUMMIT TRAVELER is charged by a USB charger, once get full charged, can continuously use for 4 hours.
- 4.SUMMIT TRAVELER adopt red 660nm and infrared 850nm=7:8, These are the 2 most efficient wavelength that are good for human health.
- 5.There are 3 different treating mode, only 660, only 850, and 660 850 half and half. You can choose the treatment mode according requirement.
- 6.SUMMIT TRAVELER use pure aluminum shell design, and adopt silent heat dissipation, so the device may get little hot when working for a while but please don't worry, this is normal.
- 7.SUMMIT TRAVELER use aluminum shell and equipped with a handle belt, suitable for one-hand use.
- 8.SUMMIT TRAVELER use slow start design. LEDs will gradually reach normal brightness after you started device, effectively avoid harm to eyes caused by sharp-flash.

Accessories

1x SUMMIT TRAVELER	1x User Manual
1x USB adapter	1x eye protector
1x travelling storage box	1x USB power cord
1x storage bag	

What are the benefits?

The benefits of red light: Red light between 600-700nm offers benefits for skin texture and tone, smoothing fine lines and wrinkles, promoting collagen production and generally rejuvenating the appearance of skin and hair.

The benefits of near-infrared light: Near-infrared light (NIR) 700nm - 1100nm. penetrates deeper into the tissues to assist with wound healing, muscle recovery, nerve injury and joint pain.

How does it work?

Red light boosts cellular energy production

The most well studied mechanism of action surrounding red and near infrared light therapy is increased mitochondrial energy production in the cells. Mitochondria are tiny organelles that are the energy power plants of all the cells in our body. Red and near infrared light therapy helps the mitochondria create more of that energy.

The specific photons found in red and near infrared light interact with a photoreceptor within our cells called cytochrome c oxidase. This interaction stimulates the mitochondria in our cells to use oxygen more efficiently, which allows the mitochondria to produce more ATP (adenosine triphosphate). Adenosine triphosphate is a complex organic chemical that provides energy to drive many processes in living cells, e.g. muscle contraction, nerve impulse propagation, and chemical synthesis. Found in all forms of life, ATP is often referred to as the "molecular unit of currency" of intracellular energy transfer. When our cells have more energy, they simply perform better and the body follows suit!

SUMMIT TRAVELER is designed as the portable LED therapy light device, it's best for target treatment use. As a adjuvant therapy device, this product cannot completely replace drugs. With SUMMIT TRAVELER, you can enjoy natural LED therapy light treatment at home. Tips here is to show how:

Before You Start:

For optimal results and safe use, follow these 3 steps before starting a treatment with the SUMMIT TRAVELER.

Step 1

Read the User Manual: Please read this instruction in full before using the SUMMIT TRAVELER.

Step 2:

Turn on the device, confirm the remaining battery capacity is enough for 1 course of treatment.

Step 3:

Set treatment mode, 660nm only, 850nm only or 850nm 660nm half and half.


Step4:

Position your devices 7,5-20 cm away from your treatment area and then enjoy your LED therapy treatment.

Starting Your Treatment:

- 1: Turn on/off Your Light.  Touch button for 5 seconds to turn on or turn off the device.
- 2:  Touch button to set the wavelength. LED screen will show the wavelength status (R- only 660nm LEDs are working, I- only 850nm LEDs are working, R I- all LEDs are working).
- 3: Three Treatment modes choose.  Touch button , LCD screen will show LED therapy treatment mode, continue to touch to choose other mode.

L1- 850nm and 660nm work together, L2-All 660nm Red Light, L3-All 850nm Infrared Light
- 4: Timing. When device start working, default 10minutes treating time. SUMMIT TRAVELER will power off automatically when countdown is over.
- 5: Battery. There are 4 indicator lights under screen to show the remaining battery capacity.

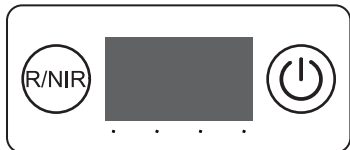Model	SUMMIT TRAVELER
LED quantity	15pcs
Irradiance	L3 mode: >150mW/cm ² at 7,5cm distance
Wavelength	660nm : 850nm= 7:8
Beam angle	30 degree
Battery capacity	8000 mAh
Dimension	L265×W160×H20mm
N.W.	0.45 KG
Warranty	2 years

This Led Light Therapy Lamp use 660nm and 850nm light wavelengths, 850nm is infrared light which is not visible by human eyes.

Maintenance

SUMMIT TRAVELER does not require any maintenance beyond basic care and regular cleaning. Cleaning the device is easy, just unplug the device and wipe the surface with a damp cloth or paper toweling. Do not use harsh cleaning chemicals or bleach because they may damage the device. Make sure the SUMMIT TRAVELER is dry before using it again. For safety, the SUMMIT TRAVELER is only intended for use indoors.

Customer Service

We are committed to your well-being by providing you with efficient, effective, and affordable light therapy devices. If at any time you are having trouble with your SUMMIT TRAVELER light therapy device, or you have any questions or concerns about the device, please do not hesitate to contact us and we will be more than happy to answer any questions you may have.

Warranty

1. Our products are warranted against defects in materials and workmanship for a period of TWO (2) YEARS from the purchase date.
2. Free return is guaranteed for a period of one MONTH.
3. Warranty does not cover the following damage(s): product been modified, improper installation or use, improper clean, electrical surges, or nature caused damage.
4. Warranty does not cover the loss of third party with direct or indirect damage(s)

To preserve the life and functionality of the SUMMIT TRAVELER, always store in a cool, dry place when possible. For proper use and safety, please be sure to read this manual and follow all instructions and do not do any of the following with SUMMIT TRAVELER:

- DO NOT let children use the SUMMIT TRAVELER without adult supervision.
- DO NOT use this device on infants.
- DO NOT let sleeping, disabled, or unconscious persons use the device.
- DO NOT exceed 20 minutes of use in one therapy session.
- DO NOT use on open, fresh wounds.
- DO NOT use in combination with liniments, salves, ointments, or balms that contain heat producing ingredients. A skin burn could be the result of doing so.
- DO NOT modify or attempt to repair this equipment. There are no parts that are serviceable by the user.

SUMMIT TRAVELER is an electrical device. To avoid electric shock and other electricity-related dangers, adhere to the following instructions:

- DO NOT wash electrical parts with water or other liquids.
- DO NOT drop your device in water or other liquids or place it where it may fall into water or other liquids. Do not expose to running water. If your device comes into contact or becomes submerged into water or other liquid, unplug it immediately from the electrical power outlet. If device is to be used in the vicinity of water, we recommend using an electrical outlet protected by Ground Fault Circuit Interrupter.
- DO NOT directly touch the device or touch the plugs or switches with wet hands.
- DO NOT apply or touch the device to, or with, wet or damp skin.
- DO NOT pull or carry the device by its cord. If the cord is damaged, don't use the device.
- DO NOT use on or near heated surfaces.
- DO NOT operate this device in areas where it could be exposed to flammable or combustible products or vapors. Explosion of fire may occur.
- DO NOT store your device in hot temperatures exceeding 140 degrees Fahrenheit.
- DO NOT use your device if it is damaged. Continuous use of a damaged device may result in electric shock or injury.
- DO NOT use extension cords unless they are designed to carry the total wattage of the device being powered.
- DO NOT unplug by pulling directly on the cord. Grasp the plug itself, not the cord, when unplugging from an outlet.
- DO NOT cover the cooling fans while operating the device.

Cautions

Recommends consulting with your healthcare provider prior to use if you believe that you are sensitive to light. Some cold and allergy medications, pain medications, and medications used to treat infection may cause light sensitivity. We also recommend consulting with your healthcare provider prior to using the SUMMIT TRAVELER if you are pregnant, have any suspicious or cancerous lesions, or have recently had a steroid injection or use topical steroids. If you experience discomfort or have a concern about the device, stop use immediately and contact your healthcare provider. Unplug from outlet when not in use to avoid tripping hazards.

Disposal

The symbol to the right, found on the SUMMIT TRAVELER or within its packaging, indicates that this device may not be disposed of with other household waste. This device requires disposal via a designated collection point for the recycling of electrical waste. Please recycle your equipment at the time of disposal to help conserve natural resources and protect human health and the environment.

Disclaimer

This product is not intended to diagnose, cure, or prevent disease. We makes no claims or representations or otherwise warranties regarding the ability of this product to cure any physical, skin or mental conditions from using this product. A qualified health professional should always be consulted prior to using this product with regards to any condition that requires medical attention

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